Hommie User Manual

Muscle Massage Gun

English [01/04]

Français[09/12]

Español[17/20]

Deutsch[05/08]



Italiano [13/16]

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Thank you for purchasing the Hommie product. Please read this [User Manual] carefully before using the product, and retain it in for future use.

■ Safety instructions and Disclaimers

1. To reduce the risk of malfunction, electric shock, injury, fire, death, and equipment damage, please read all instructions and warnings on this product in this manual before using this product.
2. For adult use only. If you are injured, it is forbidden to use this device. It can be used for recovery the muscle damage, but it needs to consult your doctor before using this device.
3. It is forbidden to use on non-dry, clean skin surface. It can be used across the clothes and gently pressed and moved on the skin for about 60 seconds.
4. Use the product only on the soft tissues of the body without causing pain or discomfort. Do not use it above the neck, any parts with rich vascular (such as neck, groin, or armpit), or any hard or boned part of the body (such as the tibia on the anterior side of the calf and the spinous process of a vertebra).

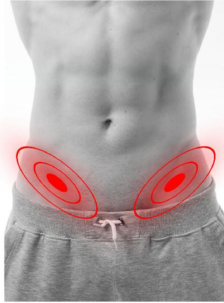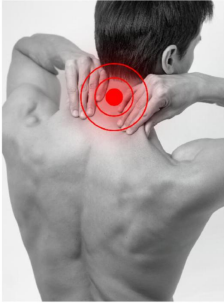5. Do not use the device in the massaging area with muscle strain or ligament injury.
6. Only use the provided best-of-breed applicator.
7. Bruises can occur regardless of restrainedly compressing or pressing. Please check the massaging area regularly and stop immediately if there is a first sign of pain or discomfort.
8. Keep your fingers, hair, or other parts of your body away from the the shaft and back of applicator as they may get caught.
9. Only use the supplied adapter to charge the device.

■ Product Overview

☑ What is Muscle Massage guns?

Muscle Massage Gun or Handheld Percussion Massager, is designed to aid in the increase of blood flow and the myofascial release process. These handheld device, use vibration and pulsation to apply concentrated pressure deep into your muscle tissue and promote recovery. This form of deep-tissue massage can work out knots and tension out of the muscles, while also enhancing your range of motion and increasing blood flow. Percussion massager, when used correctly, can help reduce your recovery time and help relieve pain, tension, and soreness.

☑ Buttons and Controls

1. Changeable applicator
2. Battery detach button
3. Speed display
4. Power display
5. Speed-
6. Speed button
7. Speed+
8. Battery LED indicator
9. Charging slot
10. Power switch

■ Operating Instructions

☑ Charging.

1. Connect the DC end of the supplied 24V adapter to the charging port of the device and plug the AC side cable into the wall outlet. **Note: Charge the device for six hours before first use.**

2. When the device is charging, the battery LED Indicator flashes red and then flashes green. After a full charge and during using, the indicator is solid green. When the battery power is low, the indicator is solid red to remind you to recharge the device.

Tips: The battery can be charged anytime. It is recommended not to completely drain the battery before recharging.

☑ Power On/Off.

Move the power switch to On position to turn on the device; Move the button to Off position to turn off.

☑ Operation.

1. Connect the the massage applicator to the device according to your needs.
2. Turn on the power switch which is at the bottom of the massager handle, and then the battery LED indicator will be solid green.
3. Press the the speed button to activate this device. You can adjust the speed level by pressing the speed +/-.

Tips: Speed 1-4: Mild Massage; Speed 5-8: Professional Massage; Speed 9-12: Relaxing Your Muscle; Speed 13-16: Decomposing Lactic Acid; Speed 17-20: Professional Mode.

■ 6 Interchangeable Applicators

		
<p>Suitable for massaging Arms, Waist, Back, Hip, Thigh, Shank and Large muscle groups.</p>	<p>Suitable for massaging Neck, Spine, Achilles.</p>	<p>Suitable for massaging the most Muscle groups, Arms, Back, Waist, Legs.</p>
		
<p>Suitable for relaxing and moulding the muscle of each part.</p>	<p>Suitable for impacting Deep-tissue, Meridians, Joints, Palm, Pelma.</p>	<p>Suitable for massaging the acupoint of inner deep-tissue, etc.</p>

☑ **Assemble and Disassemble:** When the device is off, straightly insert the desired applicator into the hole of black connector entirely. If you want to disassemble the applicator, please straightly pull the applicator out entirely.

Note: Don't rotate the applicator when assembling or disassembling the device.

■ Specifications

Speed	20 speeds; 1800-3200 percussion/minute
Input Voltage	110-240V; 30-50Hz
Power	20-30W
Amplitude	10mm
Battery	rechargeable lithium ion battery(24V/2500mAh)
Charging time	2 hours
Working time	4 hours

■ Maintenance

1. Use a slightly damp towel to clean the device and a dry with a soft cloth.
2. For storage or travel, store the device in supplied portable bag.

■ Caution and Warnings

1. Do not immerse it in water or allow water to enter the vents of the product.
2. Do not disassemble the screws or attempt to disassemble the device and modify the fascia gun in any way.
3. Do not use it continuously for over 1 hours. Using for some time, the device should be used again after 30 minutes.
4. Unplug the adapter from the device after charging or before using.
5. It is forbidden to use or charge the device with no one watching.

■ Important Reminders

1. **Don't strike your joint.** Generally, the device is suitable the soft tissues of the body. If you use it in your joint, it may damage joints.
2. **Not all parts are suitable for use.** Some parts with thin muscle and near the major artery and organs, such as neck, armpit, thorax and abdomen, etc. are not recommended to use this device.
3. **Long-time using and strong pain don't ensure the effectivity of the device.** When using, the pain should keep appropriate and the using time on the same position should be 5-10 minutes.

■ Contact Us

If any questions, please feel free to contact us: hommiehk@hotmail.com

Manufacturer: Yongkang Jiaqing Fitness Equipment Co., Ltd

Address: 2nd Floor, Chenxu Lock Industry, Kangzhuang Road, 330 National Road, Lieqiao Industrial Zone, Yongkang, Zhejiang, China

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